

Product Information

ZwickRoell Vibrophore 1000



Vibrophore 1000

CTA: 102488 154959




The ZwickRoell Vibrophore represents the latest generation of high-frequency pulsators, combining for the first time an electro-magnetically excited dynamic testing machine with a fully functional static materials testing machine. This is made possible by the 'two in one' function developed by ZwickRoell.

As well as being employed for quasi-static tensile and compression tests, Vibrophores are used in particular for determination of fatigue strength with regard to fatigue life and fatigue limit.

This allows highly efficient fatigue testing, for example to DIN 50100 (S-N curve), with tensile, compression, pulsating and alternating loads.

Other typical applications include fatigue tests on components (e.g. connecting rods, crankshafts, chains and threaded connectors) and standard metal specimens, together with fracture mechanics investigations on CT and SEB specimens,

All tests can be force, displacement or strain-controlled, while with the addition of appropriate devices testing can also be carried out under various environmental conditions (e.g. temperature, aggressive media), together with torsion and flexure tests.

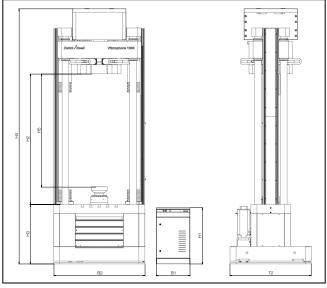

Drawing: Vibrophore 1000

Advantages and features

- fully functional static and dynamic testing machine with stiff four-column load frame
- high test-frequencies, giving short test times and high specimen throughput
- easy to operate via testXpert III and testXpert Research test programs matched exactly to the testing situations
- resonance drive has very low energy consumption (only around 2% of that of servo-hydraulic testing machines)
- easy to install no ancillary units or systems (e.g. hydraulics, coolant) needed
- maintenance-free system wear-free components
- clamping table at convenient working height for operator-friendly testing
- remote control with color display showing measurement channels, machine and test status; also used for machine set-up
- testControl II digital measurement and control electronics with 10kHz control frequency and 24-bit resolution
- continuously controlled servo motor for fast, accurate mean force control
- pulse-width modulation for high control stability and low failure liability
- easy test-frequency change simply by varying weights.
- large range of standard tools and fixtures



Product Information

ZwickRoell Vibrophore 1000

Technical data

Description	Value	
Fmax	±1000	kN
Mean load max.	±1000	kN
Force amplitude max.	±500	kN
Specimen elastic strain (oscillation displacement) max.	4 (±2)	mm
Test frequency range ¹⁾	30 285	Hz
No. of frequency steps ²⁾	8	
No. of guide columns	4	
No. of lead-screw drives	2	
Frame stiffness at 1000 mm crosshead separation ³⁾	approx. 880	kN/mm
Drive		
Static tests and for mean force control	AC drive	
min. crosshead speed	0.0001	mm/min
max. crosshead speed	250	mm/min
max. crosshead return speed	1000	mm/min
Repeat positioning accuracy on crosshead	± 2	μm
Motor holding brake	Yes	
Dynamic tests	Wear-free electro-magnetic drive	
Load-frame dimensions		
H4 - height	4538	mm
H3 - clamping-table height	1058	mm
H2 - test area height	2320	mm
H5 - test area height including load cell	2005	mm
B2 - width	1620	mm
T2 depth	1445	mm
Weight approx. ⁴⁾	17500	kg
Test area dimensions		
Test-area height without load cell	2310	mm
Test area depth	1995	mm
Free access width between columns	982 ⁵⁾	mm
Max. travel, machine frame		mm
Item No.		
Optional: Two in one for static tests	1017257	
Optional: Oscillation displacement measurement	1012273	

¹⁾ The testing frequency is determined by the stiffness and mass of the test arrangement.

⁵⁾ without bellows boot

testControl II measurement and control electronics	
Control frequency	10 kHz

²⁾ via activation / deactivation of weights

Load frame stiffness: This value is obtained via direct deformation measurement between the crossheads (moving and base crossheads) and does not take into account deformation of the drive and load cell.

⁴⁾ Weight without specimen grips or fixtures



Product Information

ZwickRoell Vibrophore 1000

testControl II measurement and control elect	ronics		
Measured-value acquisition	10 kHz, 24 bits, arithmetic	10 kHz, 24 bits, arithmetical	
Slots	5 x module bus (2 occupie	5 x module bus (2 occupied as default)	
PC interface	GigaBit Ethernet	GigaBit Ethernet	
Integrated safety concept	2-channel specification for	2-channel specification for maximum safety	
	interface for interlocked sa	interface for interlocked safety doors	
	interface for Emergency S	interface for Emergency STOP loop	
Display-equipped remote control unit	set-up or testing mode	set-up or testing mode	
	Emergency STOP button	Emergency STOP button	
	key-switch for switching b	key-switch for switching between setup and testing modes	
Dimensions, testControl II measurement and	control electronics		
H1-height without table plate	1000	mm	
B1-width	600	mm	
Approx. weight	135	kg	
Length of cable Vibrophore - testControl II	5	m	
Accessories			
Universal measurement amplifier (029443)	Choice of AC/DC supply	Choice of AC/DC supply	
	Half and full-bridge strain-	Half and full-bridge strain-gages	
	4 and 6-conductor techno	blogy	
I/O card (029448)	1 analog ± 10V input (con	1 analog ± 10V input (controllable)	
	2 analog \pm 10V outputs	2 analog ± 10V outputs	
	4 digital inputs, 24V		
	3 digital outputs, 24V	3 digital outputs, 24V 1 relay output, potential-free	
	1 relay output, potential-fr		

Installation conditions		
Operating temperature	+10 +35	°C
Storage temperature	-25 + 50	°C
Humidity (non-condensing)	20 90	%
Electrical supply		
power supply voltage 3 Ph/N/PE	400	V
power frequency	50 / 60	Hz
back-up fuse	16	А
Noise level at 1 meter distance ¹⁾	80 110	dB(A)

¹⁾ Depends on test arrangement, test load and testing frequency